



**Megazyme,**  
Bray Business Park,  
Bray, Co. Wicklow,  
A98 YV29,  
Ireland.  
Tel: + 353 1 286 1220  
Fax: + 353 1 286 1264

## **Validation Report: Ammonia Assay Kit (Rapid) (cat. no. K-AMIAR)**

### **1. Scope**

Megazyme's Ammonia Assay Kit (Rapid) (K-AMIAR) is an enzymatic method used for the rapid measurement and analysis of ammonia in foodstuffs, beverages and other materials. This method was developed in-house and measures ammonia in g/L. Methods based on this principle have been accepted by MEBAK.

### **2. Planning**

The purpose of this report is to verify and validate the current method as detailed by Ammonia Assay Kit (Rapid) (K-AMIAR).

### **3. Performance characteristics**

The selectivity, working range, limit of detection, limit of quantification, trueness (*bias*) and precision of this kit is detailed in this report.

#### **3.1. Selectivity**

This assay is specific for ammonia. Interfering substances in the sample being analysed can be identified by including an internal standard. Quantitative recovery of this standard would be expected. Losses in sample handling and extraction are identified by performing recovery experiments, i.e. by adding ammonia to the sample in the initial extraction steps.

#### **3.2. Working Range**


Assay follows the Ammonia Assay Kit (Rapid) (K-AMIAR) standard procedure. 0.1 mL of Ammonia standard was used as sample, with a range of concentrations (0.002-0.07 g/L ammonia) which corresponds to 0.2-7.0 µg of ammonia per cuvette. Absorbance A<sub>2</sub> was read after 5 min, at 340 nm and at 25°C as recommended in the procedure.



**Megazyme,**

Bray Business Park,  
Bray, Co. Wicklow,  
A98 YV29,  
Ireland.  
Tel: + 353 1 286 1220  
Fax: + 353 1 286 1264

Ammonia Concentration [ $\mu\text{g}/\text{cuvette}$ ]	$\Delta A_{340\text{nm}}$
0	0.000
0.2	0.026
1	0.132
2	0.276
3	0.413
4	0.551
5	0.696
7	0.977

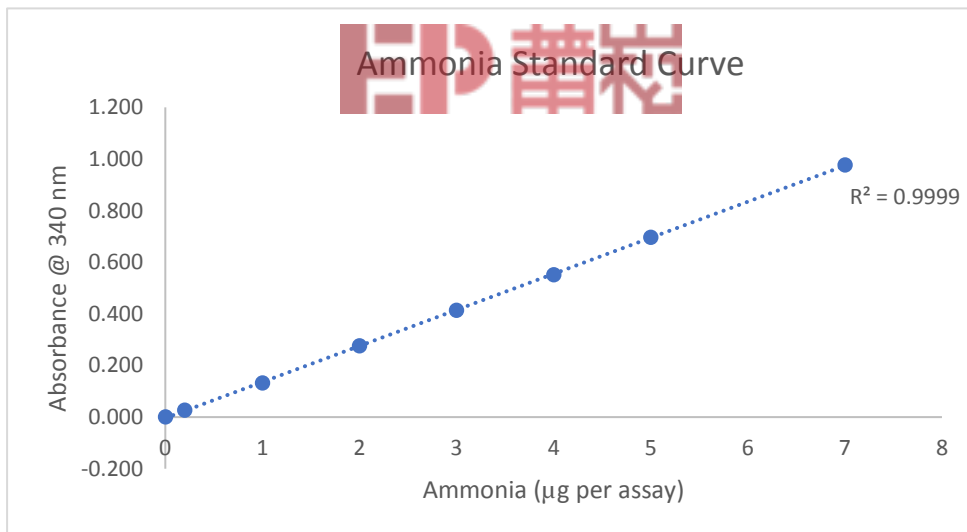



**Megazyme,**

Bray Business Park,  
Bray, Co. Wicklow,  
A98 YV29,  
Ireland.  
Tel: + 353 1 286 1220  
Fax: + 353 1 286 1264

### 3.3. LOD and LOQ

The **instrument limit of detection**, as per kit booklet, is 0.071 mg/L, which is derived from an absorbance difference of 0.020 with the maximum sample volume of 2.00 mL.

The **calculated limit of detection (LOD)** and the **calculated limit of quantification (LOQ)** for this report purpose is based on the analysis of samples that have been taken through the whole Ammonia Assay Kit (Rapid) (K-AMIAR) measurement procedure.

- The LOD is the lowest concentration of the analyte that can be detected by the method. LOD is calculated as  $3 \times s'0$ ; where  $s'0$  is the standard deviation of a number of samples A1 reading.
- The LOQ is the lowest level at which the kit's performance is acceptably repeatable. LOQ is calculated as  $kQ \times s'0$ ; where  $s'0$  is the standard deviation of a number of samples A1 reading. The IUPAC default value for  $kQ$  is 10
- For Ammonia Assay Kit (Rapid) (K-AMIAR)

**LOD – For 2.0 mL of sample (maximum volume)**

Ammonia = 0.050 mg/L

**LOQ – For 2.0 mL of sample (maximum volume)**

Ammonia = 0.150 mg/L

\* **Note:** The above detection limits are for samples as used in the assay, after sample preparations if required (e.g. deproteinisation). The dilution used in pre-treatment must be accounted for while establishing the detection limits for specific samples.



**Megazyme,**

Bray Business Park,  
 Bray, Co. Wicklow,  
 A98 YV29,  
 Ireland.  
 Tel: + 353 1 286 1220  
 Fax: + 353 1 286 1264

**3.4. Trueness (Bias)**

Comparison of the mean of the results (x) achieved with Ammonia Assay Kit (Rapid) (K-AMIAR) method with a suitable reference value (x ref). For this report, Relative Bias is calculated in per cent as:  $b(\%) = \frac{x - x_{ref}}{x_{ref}} \times 100$ . The reference material for this purpose is ammonia supplied with the Ammonia Assay Kit (Rapid) (K-AMIAR) at 0.04 g/L.

**Relative Bias b(%)**

	n	Ref Material (g/L)	Mean (g/L)	b(%)
<b>Ammonia</b>	22	0.04	0.0396	-1.06

**3.5. Precision**

This report details the reproducibility of the Ammonia Assay Kit (Rapid) (K-AMIAR), it is a measure of the variability in results, on different days and by different analysts, over an extended period of time. For the purpose of this report different lot numbers of the kit standard is used as the reference material.

**Reproducibility**

	n	Ref Material (g/L)	Mean (g/L)	Standard Deviation	%CV
<b>Ammonia</b>	22	0.04	0.0396	0.0005	1.29

Repeatability of this kit can be assessed using wine samples. This is a measure of the variability in results by a single analyst, using real samples, using the same equipment and over a short period of time. The use of wine samples shows one of the many applications of this kit.

**Repeatability**

	n	Mean (g/L)	Standard Deviation	%CV
<b>White Wine</b>	12	0.0091	0.0002	2.09
<b>Red Wine</b>	12	0.0220	0.00035	1.59



**Megazyme,**

Bray Business Park,  
Bray, Co. Wicklow,  
A98 YV29,  
Ireland.  
Tel: + 353 1 286 1220  
Fax: + 353 1 286 1264

**4. Conclusion**

The method outlined in this document is a robust, quick and easy method for the measurement of ammonia in various matrices. It has been used for many years and is fully automatable for high throughput analysis of samples. Data presented in this report verifies and validates that this method is fit for the purpose intended, which is summarised below.

<b>Validation Summary</b>	<b>Ammonia</b>
<b>Working range (<math>\mu\text{g}</math> in cuvette)</b>	0.2-7.0
<b>LOD (mg/L)</b>	0.05
<b>LOQ (mg/L)</b>	0.15
<b>Relative Bias <i>b</i> (%)</b>	-1.06
<b>Reproducibility (%CV using kit standard)</b>	1.29
<b>Repeatability (%CV using white wine)</b>	2.09
<b>Repeatability (%CV using red wine)</b>	1.59

